



Passport™ BTE Guide

Moxi™ 13 / Moda™ 13



Table of Contents

Your Passport™ Hearing Instruments	1
Overview of your Passport™ Hearing Instrument	2
Features of your Passport™ Hearing Instrument	4
Putting your Hearing Instruments on your Ears	5
Turning your Hearing Instruments On and Off	6
Battery Information.....	7
Operating Instructions	9
Caring for your Hearing Instruments.....	16
Cleaning your Hearing Instruments	17
Signature Features of your Passport™ Hearing Instruments.....	20
Assistive Listening Devices	21
Warnings.....	22
Troubleshooting Guide	25
Warning to Hearing Instrument Dispensers	28

Your Passport™ Hearing Instruments

Hearing Healthcare Professional: _____

Telephone: _____

Model: _____

Serial Number: _____

Replacement Batteries: Size 13

Warranty: _____

Program 1 is the Automatic Program

Program 2 is the manual program for: _____

Program 3 is the manual program for: _____


Program 4 is the manual program for: _____

Date of Purchase: _____


Overview of your Passport™ Hearing Instrument

The diagrams below identify some of the components on your hearing instruments. Have your hearing healthcare professional place a check beside the diagram that best describes your hearing instrument style and attachment.

Moda™13 Attachment Style (check one)

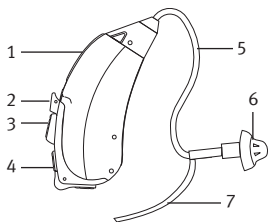

Moda 13 with slim tube and dome

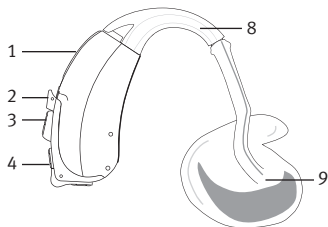


Moda 13 with earmold


Moda 13 Legend

- 1 Microphone and Microphone Shield
- 2 Battery Door/ On & Off Switch
- 3 Program Button or Volume Control (depending on your customized fitting)
- 4 Toggle for volume control or smartFocus™
- 5 Tube
- 6 Dome
- 7 Retention Piece
- 8 Earhook
- 9 Earmold

Moxi™13 Attachment Style (check one)

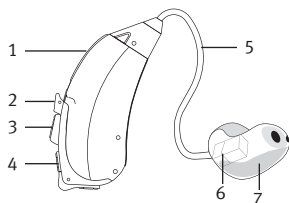

Moxi 13 with Shell

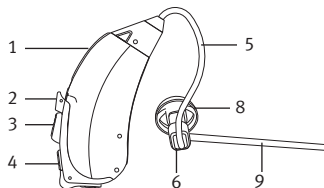

Moxi 13 with Dome

Moxi 13 Legend

- 1 Microphone and Microphone Shield
- 2 Battery Door/ On & Off Switch
- 3 Program Button or Volume Control (depending on your customized fitting)
- 4 Toggle for volume control or smartFocus™
- 5 Tube
- 6 Speaker
- 7 Shell
- 8 Dome
- 9 Retention Piece

My hearing instrument has (check all that apply)

- | | |
|--|---|
| <input type="checkbox"/> smartFocus™ Toggle | <input type="checkbox"/> Push Button Program Button |
| <input type="checkbox"/> Volume Control Toggle | <input type="checkbox"/> Push Button Volume Control |
| <input type="checkbox"/> duoLink | |

Optional Accessories (check all that apply)

- | | |
|--|-----------------------------------|
| <input type="checkbox"/> Smart Control | <input type="checkbox"/> uDirect™ |
|--|-----------------------------------|

See Smart Control and uDirect User Guides for more information.

Features of your Passport™ Hearing Instrument

Passport premium hearing instruments come equipped with the following features:

- automatic listening program
- 3 additional manual programs
- SmartFocus with personal adjustments for Clarity and Comfort settings
- Wireless capability with duoLink
- Optional wireless connectivity to other devices


It may take some time for you to fully adjust to your new hearing instruments. In the beginning, do not feel as though you need to wear them for longer than is comfortable. Depending on your previous experience with hearing instruments you can increase wearing time gradually until you have fully adjusted to them. Eventually, you should wear your passport hearing instruments all day every day.

Please consult your hearing healthcare professional regarding any adjustments you may require to ensure that you are happy with your hearing instrument. For instance, pitch and loudness of the beeps can be adjusted or turned off entirely. Also note that this is a general guide and your hearing healthcare professional may have enabled or disengaged certain features depending on your requirements.


Putting your Hearing Instruments on your Ears

Your Passport hearing instruments may be color-coded red for your right ear and blue for your left ear. This color indicator is located inside the battery door with a small colored dot.

Passport™ with Domes

1. Hold the tube where it attaches to the dome and gently push the dome into your ear canal. The tube should lie flush against your head and not stick out.
2. Place Passport over the top of your ear.
3. Place the retention piece in your ear so it rests at the bottom of the opening of your ear canal.

Passport™ with Earmolds/Shells

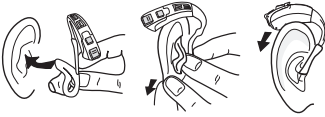1. Carefully insert the earmold/shell into your ear. The earmold/shell should fit into your ear snugly and comfortably. The tube should lie flush against your heard and not stick out.
2. Place Passport over the top of your ear.

Turning your Hearing Instruments On and Off


Your hearing instruments have a three-position battery door that acts as an off/on switch and that allows access to the battery compartment.

ON:

To turn the hearing instruments on, close the battery door fully. *Note: It may take 5 seconds before the hearing device turns on. Your hearing healthcare provider can additionally increase the start up delay if required.*

OFF:

To turn the hearing instrument off, partially open the battery door. This position also allows excess moisture to vent away from the battery compartment.

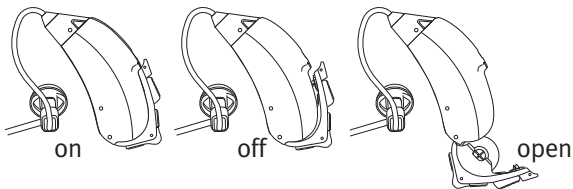

Note: When turning your hearing instrument on and off while it is on the ear, grasp the top and bottom of the instrument with your index finger and thumb. Use the index finger of your opposite hand to open and close the battery door.

Battery Information


To replace the hearing instrument battery, fully open the battery door for access to the battery compartment.

Low Battery Warning

When you hear two long beeps, your hearing instruments are warning you that their batteries are low. You will hear the warning approximately every 30 minutes until you change the batteries or the batteries die. After the first warning, you may experience some reduction in sound quality. This is normal and can be remedied by inserting fresh batteries in the hearing instruments.

Replacing the Battery

1. Gently swing out the battery door with your fingernail, or grasp the top and bottom of the instrument with your index finger and thumb. You can use your index finger to open and close the battery door if this is easier.
2. Grasp the battery with your thumb and index finger and remove.

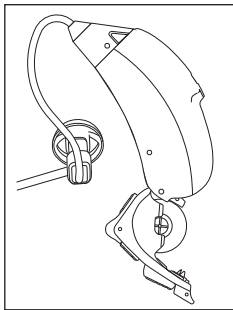

3. Insert the new battery into the battery compartment matching it with the battery symbol inside the battery door. This will ensure that the battery door closes properly.

Note: If the battery is inserted incorrectly, the door will not close.

4. Close the battery door.

Caring for Batteries

- Always discard batteries in a safe and environmentally friendly way.
- To prolong battery life, remember to turn your hearing devices off when not in use.
- Remove the batteries and keep the battery door open while hearing devices are not in use. This will allow internal moisture to evaporate.

Operating Instructions

Push Button








Your hearing instrument is equipped with a push button, built into the back of the battery door, which has been set for switching between programs or adjusting your volume. Note that your push button may be disabled by your hearing healthcare professional and, therefore, could be inactive.

Program Control

If you have an active push button program control, each time you push the button, you will move to a new program. Also if duoLink is enabled, the hearing instrument is set so that by pressing the push button on one hearing instrument, both hearing instruments will change programs.

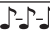Your hearing instruments beep to indicate which

Program 1 (e.g. Automatic Program)		1 beep
Program 2 (e.g. Group/Party Noise)		2 beeps
Program 3 (e.g. Easy-t/Telephone)		3 beeps
Program 4 (e.g. Music)		4 beeps
Easy-t/Telephone		short melody
duoLink (if enabled)		Opposite ear beeps the same as side being adjusted

program you are in (i.e. one beep for program 1, two beeps for program 2, etc.)




Volume Control

If your push button has been configured as a volume control, you can adjust the volume level by pushing the button.

If duoLink is enabled, the hearing instrument can also be set so that by changing the volume on one hearing instrument, both will be adjusted.

Your Passport hearing instruments will learn your volume control adjustments while the hearing instruments are set to the automatic program, if this feature has been enabled by your hearing healthcare professional. Over time, the hearing instruments will adjust the default volume settings in the automatic program to the levels that you typically use.

Your hearing healthcare professional can help you identify your different push button volume control settings.

Volume Setting	Beeps
Ideal volume level	 1 beep
Increased loudness	 1 beep + 1 high-pitched beep
Decreased loudness	 1 beep + 1 low-pitched beep

Your Passport hearing instruments may also come with an optional remote control which allows you to switch between different listening programs, and adjust your smartFocus. Please refer to your Smart Control User Guide for more information.

Toggle



◦ The toggle on your hearing instruments can be set up as either a volume control, or a smartFocus control. Note that your toggle may be disabled by your hearing healthcare professional and, therefore, could be inactive.




Volume Control

If your toggle has been configured as a volume control, you can increase the volume level by pushing in the toggle portion on the back of the hearing instrument. You can decrease the volume level by pushing the toggle portion located on the bottom of the hearing instrument.

If duoLink is enabled, the hearing instrument can also be set so that by changing the volume on one hearing instrument, both will be adjusted.

Your Passport hearing instruments will learn your volume control adjustments while the hearing instruments are set to the automatic program, if this feature has been enabled by your hearing healthcare

professional. Over time, the hearing instruments will adjust the default volume settings in the automatic program to the levels that you typically use.

Volume Setting	Beeps
Ideal volume level	 1 beep
Increased loudness	 1 beep + 1 high-pitched beep
Decreased loudness	 1 beep + 1 low-pitched beep




Your hearing healthcare professional can help you identify your different volume control settings.

SmartFocus™

If the toggle has been configured as a smartFocus, push in the toggle portion on the back of the hearing instrument to improve the clarity of sounds in front of you, such as speech. For greater overall listening comfort in noisy listening situations, push the portion of the toggle located at the bottom of the hearing instrument. Each step will change the hearing instrument performance. Sometimes multiple steps are required to achieve the desired listening results. After each step wait for about 4-5 seconds to allow the hearing instrument to adjust to the new setting.

If duoLink is enabled, the hearing instrument can also be set so that by changing the volume on one hearing instrument, both will be adjusted.

As you change the smartFocus level, your hearing instruments will beep. Please refer to the chart below to learn more about the different comfort-clarity beeps.

SmartFocus Setting	Beep
Maximum sound clarity level	 2 beeps
Maximum listening comfort level	 2 beeps
Ideal comfort and clarity level	 1 beep

Using the Telephone

As telephones do not all work the same, you may experience different results from different phones. You can use many phones by simply holding the receiver to your ear, without the need to change to a dedicated telephone program. Remember to move the handset slightly up or back to find the position that sounds best to you.

Depending on the phone type you use, your hearing healthcare professional might have selected a dedicated phone program on your hearing instruments. Therefore your hearing instrument may have an automatic telephone program (easy-t feature), which automatically switches you to a dedicated telephone program when your hearing instrument is close to a telephone receiver. When the receiver is moved away from the hearing instrument it will

automatically return to the previous listening program. If the hearing instrument does not switch to the telephone program automatically when the telephone receiver is placed in proximity, the magnet for easy-t hearing instruments should be attached to the telephone receiver. The magnet is designed to strengthen the magnetic field at the ear piece of hearing instrument compatible telephones.

If your Passport hearing instruments are configured to work with the wireless uDirect accessory, you might benefit from a Bluetooth hands free option. Refer to the uDirect User Guide for more information.

If duoLink is enabled for the automatic phone program, both hearing instruments are synchronized and the listening preference for the opposite ear will also change when an automatic phone program is activated.

Consult your hearing healthcare professional if you experience difficulties using your hearing instruments on the phone.

To affix the optional easy-t magnet:

1. Clean the telephone receiver.
2. Hold the magnet near the “listening end” of your telephone receiver and release it (Figure 1). The magnet will flip to the appropriate side and seek

the optimal position on the telephone receiver.

3. Place the double-sided tape in this optimal position on the telephone receiver (Figure 2) and attach the magnet to the tape (Figure 3).

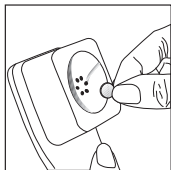


Figure 1

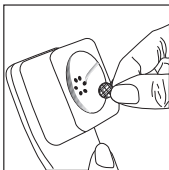


Figure 2



Figure 3

Caring for your Hearing Instruments

- Open the battery door when not in use.
- Always remove your hearing instruments when using hair care products. The hearing instruments can become clogged and cease to function properly.
- Do not wear your hearing devices in the bath or shower or immerse them in water.
- If your hearing instruments do become wet, do not attempt to dry them in an oven or microwave. Do not adjust any controls. Open the battery doors immediately, and allow your hearing instruments to dry naturally for 24 hours.
- Protect your hearing devices from excessive heat (hair dryer, vehicle glove box or dashboard).
- Ensure you do not twist or squeeze the tube when your hearing instruments are placed in their case.
- Regular use of a dehumidifier, such as a Dri-Aid kit, can help prevent corrosion.
- Do not drop your hearing devices or knock them against hard surfaces.

Cleaning your Hearing Instruments


Ear wax is natural and common. Ensuring your hearing devices, earmolds and domes are free of ear wax is an important step in your daily cleaning and maintenance routine.

- Never use alcohol to clean your hearing devices, earmolds or domes.
- Do not use sharp tools to dislodge ear wax. Sticking household items into your hearing devices or earmolds can seriously damage them.

Passport Hearing Instruments with Earmolds

Earmolds send amplified sound from the hearing devices into the ears.

They must fit into your ears snugly and comfortably. If amplified sound leaks out of your ear, you may hear whistling.

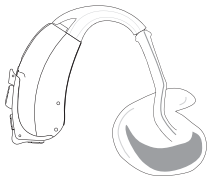


Always ensure earmolds are clean and free of earwax. Do not use alcohol to clean your earmolds. If the earmolds become plugged, clear the opening with a wax loop or pipe cleaner. If your physician prescribes eardrops, clean any moisture that may get into the earmolds or tubing to prevent plugging.

If the earmolds require further cleaning, disconnect the plastic tube from the hook of the hearing devices.

Wash only the earmolds in warm water with a mild soap. Rinse them with cool water and allow them to dry overnight. Make sure the earmold tubes are dry before reconnecting them to each hook on your hearing devices.


Passport Hearing Instruments with Slim Tubes and Domes

You should have the slim tubes and domes replaced by your hearing healthcare professional approximately every three to six months or when they become stiff or brittle. Clean the domes daily with a damp cloth. You should also clean the slim tubes periodically, with the cleaning pin provided, when you begin to notice debris in and around the tubes.

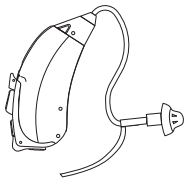

1. Hold the slim tube in one hand and the hearing instrument in the other hand.
2. Gently turn the hearing instrument until it detaches from the slim tube.
3. Use a damp cloth to clean the outside of the slim tube and dome.
4. Remove the dome from the slim tube before cleaning the slim tube.

- Using the black cleaning pin provided in the kit, gently insert the cleaning pin where the slim tube attaches to the hearing instrument and push the pin all the way through the tube.

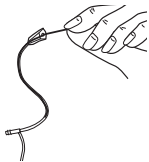


Note: The slim tubes and domes should never be rinsed or submerged in water as water drops may become lodged in the tube, block sound or damage the electrical components of the hearing instrument.

- Once the slim tube has been cleaned, reattach it by gently turning the hearing instrument onto the slim tube.

Passport Hearing Instruments with Shells and Domes

- Clean the domes and shells on the outside daily with a damp cloth. Avoid getting any water in and around the speaker units and shells.

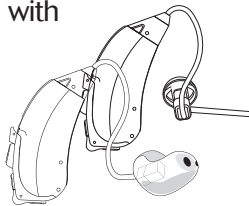

Note: The tubes, speakers, domes or shells should never be rinsed or submerged in water as water drops may block sound or damage the electrical components of the hearing instrument.

- If your shells require further cleaning, the speaker waxguard may be plugged and require replacing. See your hearing healthcare professional.

Signature Features of your Passport™ Hearing Instruments

SmartFocus™

With a Passport hearing instrument, you can adjust the focus on hearing speech or background noise in your listening environment via the smartFocus control. This control provides an additional fine tuning adjustment beyond a traditional volume control. You can adjust your smartFocus control to focus on speech (clarity) or soften background noise (comfort) using the lever on your hearing instruments or dial on your Smart Control.

As you adjust smartFocus, your Passport hearing instruments will learn your smartFocus clarity and comfort settings while the hearing instruments are set to the automatic program. Over time, the hearing instruments will adjust the pre-programmed smartFocus settings to your preferred levels.

The optional Smart Control hand-held remote control allows you to adjust your smartFocus clarity and comfort levels.

DuoLink

With duoLink, your hearing instruments wirelessly communicate with each other. As you make a program or volume change on one hearing instrument, the

other hearing instrument automatically changes as well. An additional benefit occurs while using the phone since as one hearing instrument detects the phone, the other hearing instrument can be programmed to mute sound to eliminate distractions.

uDirect™

Your Passport hearing devices may come with an optional uDirect™, a device that is worn around your neck, and provides secure connectivity between your hearing instruments and Bluetooth-enabled devices (i.e. cell phones), wired audio input jacks, or FM transmitted signals. uDirect provides easy hands-free access to clear stereo, or audio signals automatically activated when a signal is detected.

Assistive Listening Devices

Listening in Public Places

Passport's telecoil option can also help you listen in public places equipped with telecoil compatible assistive listening devices such as a loop system. When you see this symbol, it means that there is loop system installed; this loop system is compatible with your hearing instrument. Please contact your hearing healthcare provider for more information on loop systems.



Warnings

- Hearing devices should only be used as directed by your physician or hearing healthcare professional.
- Hearing devices will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions.
- Do not use your hearing devices in explosion hazard areas.
- Allergic reactions to hearing devices are unlikely. However, if you experience itching, redness, soreness, inflammation or a burning sensation in or around your ears, inform your hearing healthcare professional and contact your physician.
- In the unlikely case that any parts remain in the ear canal after the removal of the hearing instrument, contact a physician immediately.
- Remove your hearing devices for CT and MRI scans or for other electromagnetic procedures.
- Special care should be exercised in wearing hearing devices when maximum sound pressure levels exceed 132 decibels. There may be a risk of impairing your remaining hearing. Speak with your hearing healthcare professional to ensure the maximum output of your hearing devices is suitable for your particular hearing loss.

Magnet Warnings

- Be sure the magnet is securely affixed to the telephone.
- Keep loose magnets out of reach of children and pets.
- If the magnet falls into your ear, contact your hearing healthcare professional.
- If the magnet is swallowed, contact your physician immediately.
- The magnet may affect some medical devices or electronic systems. Always keep the magnet (or the telephone equipped with the magnet) at least 30 cm (12”) away from pacemakers, credit cards, floppy disks or other magnetically sensitive devices.
- Too high distortion during dialing or phoning may mean that the phone handset is stressed by the magnet. To avoid any damage, please move the magnet to another place on the telephone receiver.

Battery Warnings

- Never leave hearing devices or batteries where small children and pets can reach them.
- Never put hearing devices or batteries in your mouth. If a hearing device or battery is swallowed, call a physician immediately.

Warning to hearing healthcare professional:

- Domes should never be fitted on patients with perforated eardrums, exposed middle ear cavities, or surgically altered ear canals. In the case of such a condition, we recommend to use a customized ear mold.

Precautions

- The use of hearing devices is only part of hearing rehabilitation; auditory training and lip reading instruction may be required as well.
- In most cases, infrequent use of hearing devices does not provide full benefit. Once you have become accustomed to your hearing devices, wear your hearing devices everyday all day.
- Your hearing devices use the most modern components to provide the best possible sound quality in every listening situation. However, communication devices such as digital cell phones can create interference (a buzzing sound) in hearing devices. If you experience interference from a cell phone being used close by, you can minimize this interference in a number of ways. Switch your hearing devices to another program, turn your head in a different direction or locate the cell phone and move away from it.

Labeling

The serial number and year of manufacture are located inside the battery door.

Troubleshooting Guide

CAUSE	POSSIBLE REMEDY
No sound	
<ul style="list-style-type: none">• Not turned on• Low/dead battery• Poor battery contact• Battery upside down• Moda 13: Earmolds/slim tubes/domes blocked with ear wax• Moxi 13: Shells/domes blocked with ear wax• Plugged microphone protector	<ul style="list-style-type: none">• Turn on• Replace battery• Consult your hearing healthcare professional• Insert battery plus (+) side up• Clean earmolds or slim tubes and domes. See “Cleaning your Hearing Instrument”. Consult your hearing healthcare professional.• Consult your hearing healthcare professional
Not loud enough	
<ul style="list-style-type: none">• Low volume• Low battery• Moda 13: Earmolds/slim tubes/domes not inserted properly• Moxi 13: Shells/domes not inserted properly• Change in hearing• Moda 13: Earmolds/slim tubes/domes blocked with ear wax• Moxi 13: Shells/domes blocked with ear wax	<ul style="list-style-type: none">• Turn up volume; see hearing healthcare professional for models without a manual volume control or if problem persists.• Replace battery• See “Putting Your Hearing Instruments on Your Ears”. Reinsert carefully.• Consult your hearing healthcare professional• Clean shells and domes. See “Cleaning your Hearing Instrument”. Consult your hearing healthcare professional.

CAUSE	POSSIBLE REMEDY
<ul style="list-style-type: none"> • Plugged microphone shield 	<ul style="list-style-type: none"> • Consult your hearing healthcare professional
Intermittent	
<ul style="list-style-type: none"> • Low battery • Dirty battery contact 	<ul style="list-style-type: none"> • Replace battery • Consult your hearing healthcare professional
Two long beeps	
<ul style="list-style-type: none"> • Low battery 	<ul style="list-style-type: none"> • Replace battery
Whistling	
<ul style="list-style-type: none"> • Moda 13: Earmolds/slim tubes/domes not inserted properly • Moxi 13: Shells/domes not inserted properly • Hand/clothing near ear • Poorly fitting shells/domes 	<ul style="list-style-type: none"> • Remove and reinsert carefully • Remove hand/clothing from ear • Consult your hearing healthcare professional
Not clear, distorted	
<ul style="list-style-type: none"> • Moda 13: Poorly fitting earmolds/slim tubes/domes • Moxi 13: Poorly fitting shells/domes 	<ul style="list-style-type: none"> • Consult your hearing healthcare professional

CAUSE	POSSIBLE REMEDY
<ul style="list-style-type: none"> • Moda 13: Earmolds/slim tubes/domes blocked with ear wax • Moxi 13: Shells/domes blocked with ear wax • Low battery • Plugged microphone protector 	<ul style="list-style-type: none"> • Clean shells and domes. See “Cleaning Your Hearing Instrument”. Consult your hearing healthcare professional. • Replace battery • Consult your hearing healthcare professional

Shells/domes falling out of ear

<ul style="list-style-type: none"> • Poorly fitting earmolds/slim tubes/shells/domes • Earmolds/shells/domes not inserted properly 	<ul style="list-style-type: none"> • Consult your hearing healthcare professional • See “Putting Your Hearing Instruments on Your Ears”. Reinsert carefully.
--	--

Weak on the telephone

<ul style="list-style-type: none"> • Telephone not positioned properly • Hearing instrument requires adjustment 	<ul style="list-style-type: none"> • Move telephone receiver around ear for clearer signal. See “Using the Telephone”. • Consult your hearing healthcare professional
---	---

For any problems not listed in the guide, contact your hearing healthcare professional. If you do not have a hearing healthcare professional, please contact the nearest office listed on the back page of this booklet.

Warning to Hearing Instrument Dispensers

A hearing instrument dispenser should advise a prospective hearing instrument user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing instrument if the hearing instrument dispenser determines through inquiry, actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions: (i) Visible congenital or traumatic deformity of the ear. (ii) History of active drainage from the ear within the previous 90 days. (iii) History of sudden or rapidly progressive hearing loss within the previous 90 days. (iv) Acute or chronic dizziness. (v) Unilateral hearing loss of sudden or recent onset within the previous 90 days. (vi) Audiometric air-bone gap equal to or greater than 15 decibels at 500 hertz (Hz), 1,000 Hz, and 2,000 Hz. (vii) Visible evidence of significant cerumen accumulation or a foreign body in the ear canal. (viii) Pain or discomfort in the ear. Special care should be exercised in selecting and fitting a hearing instrument whose maximum sound pressure level exceeds 132 decibels because there may be risk of impairing the remaining hearing of the hearing instrument user. [This provision is required only for those hearing instruments with a maximum sound

pressure capability greater than 132 decibels (dB).]

Important Notice for Prospective Hearing Instrument Users

Good health practice requires that a person with a hearing loss have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing a hearing instrument.

Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists. The purpose of medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before the hearing instrument is purchased. Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for a hearing instrument. The physician will refer you to an audiologist or a hearing instrument dispenser, as appropriate, for a hearing instrument evaluation. The audiologist or hearing instrument dispenser will conduct a hearing instrument evaluation to assess your ability to hear with and without a hearing instrument. The hearing instrument evaluation will enable the audiologist or dispenser to select and fit a

hearing instrument to your individual needs. If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or purchase-option program. Many hearing instrument dispensers now offer programs that permit you to wear a hearing instrument for a period of time for a nominal fee after which you may decide if you want to purchase the hearing instrument. Federal law restricts the sale of hearing instruments to those individuals who have obtained a medical evaluation from a licensed physician. Federal law permits a fully informed adult to sign a waiver statement declining the medical evaluation for religious or personal beliefs that preclude consultation with a physician. The exercise of such a waiver is not in your best health interest and its use is strongly discouraged.

Children With Hearing Loss

In addition to seeing a physician for a medical evaluation, a child with a hearing loss should be directed to an audiologist for evaluation and rehabilitation since hearing loss may cause problems in language development and the educational and social growth of a child. An audiologist is qualified by training and experience to assist in the evaluation and rehabilitation of a child with a hearing loss.

Unitron Distributors

CORPORATE OFFICE/INTERNATIONAL

20 Beasley Drive, P.O. Box 9017, Kitchener, ON N2G 4X1 Canada

AUSTRALIA

Level 2, Norwest Quay, 21 Solent
Circuit, Baulkham Hills, NSW,
NSW 2153

BELGIUM

Baron de Vironlaan, 60
b-1700 Dilbeek

CANADA

20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1

CHINA

No. 200 Suhong Road, Export
Processing Zone 4D, Suzhou
Industrial Park 215021

DENMARK

Nitivej 10, DK-2000 Frederiksberg

EUROPEAN REPRESENTATIVE

Daimlerstrasse 22, 70736
Fellbach-Oeffingen, Germany

FRANCE

5, rue Maryse Bastié - BP 15 69671
Bron Cedex, France

NETHERLANDS

Ir. D.S. Tuijnmanweg 10,
4131 PN VIANEN

NEW ZEALAND

10/215 Rosedale Rd., M277 Private
Bag, 300987, Albany, Auckland

NORWAY

Brugata 14 0186
Oslo, Norway

SOUTH AFRICA

First Floor – Selborne House,
Fourways Golf Park, Roos Street,
Fourways, Johannesburg

SPAIN

Avenida de la Industria nº13/15,
Polígono de Canastell,
03690 Sant Vicent del Raspeig.
Alicante, España

SWEDEN

Förmansvägen 2, 4 tr
SE-11743, Stockholm, Sweden

UNITED KINGDOM

St. George House, Cygnet Court,
Centre Park, Warrington, Cheshire
WA1 1PD

U.S.A.

Suite A, 2300 Berkshire Lane North,
Plymouth, MN 55441

Manufacturer

UNITRON HEARING LTD.

20 Beasley Drive, P.O. Box 9017, Kitchener, ON N2G 4X1 Canada

CE
0543



 **unitron™**

DISTRIBUTOR

www.unitron.com

08-064 029-5735-02

